



# MEDIDENTA SURGICAL SERIES

The Air-Free and Foresight Series



## REAR EXHAUST

All air is vented through pilot holes located on the back-end. This helps provide an air free oral cavity along with great turbine ventilation for greater power.

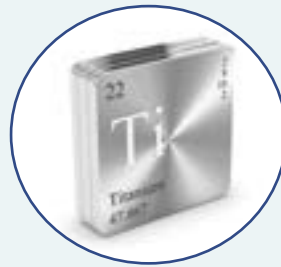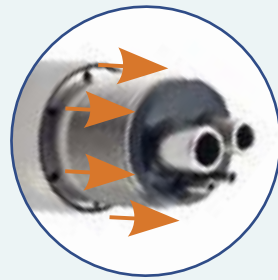

## TITANIUM FINISHED STAINLESS STEEL

The titanium finished shells are capable of withstanding the high demand of your office, preventing corrosion and scratches for long lasting operation.



## AIR FREE HEAD

Put safety first. Unlike most conventional surgical highsPEEDS, no air vents out of the head of the Air-Free handpiece. This eliminates the risk of subcutaneous air emphysema for a safe, risk-free surgical procedure.

# AIR-FREE™ 90S

## THE FIRST EVER 90° SURGICAL

The first 90-degree surgical highsPEED to accommodate the many surgical procedures that do not require a 45-degree angle, for a comfort and efficiency.



- 90° Angle
- Ideal for perio and osseous surgical procedures

**No chance for any air/oil contamination at surgery site.**

# AIR-FREE™ 45

## FOR SAFER THIRD MOLAR EXTRACTIONS



- 45° Angle
- Ideal for mesioangular impactions and third molar extractions

## AIR FREE 90S

\$549.99 | 1 YEAR WARRANTY | ITEM# AF90SG

### Air-Free Series

- 20 Watts of Power
- Pure Water Jetstream Cooling
- Ceramic Bearings
- Fixed 4-Hole
- Made in USA

### 90 Degrees for Ergonomic Comfort



## AIR FREE 45

\$599.99 | 1 YEAR WARRANTY | ITEM# AF45

### Air-Free Series

- 20 Watts of Power
- Pure Water Jetstream Cooling
- Ceramic Bearings
- Fixed 4-Hole
- Made in USA

### No chance for any air/oil contamination at surgery site.





# SELF GENERATING LED

## What is a Self Generating LED?

Self Generating LED (SGL) is a system inside the handpiece that produces LED optics without the need for optics set up in your chairline.



\$799.99 | 1 YEAR WARRANTY | ITEM# FS-45S


## FORESIGHT SGL SERIES

Traditional fiber optics are now a thing of the past. The Foresight SGL (Self Generating LED) Series offers bright, direct illumination to your procedure site. The Foresight 45 (surgical) works exceptionally well in illuminating the third molar root tips. The built-in micro generator activates the LED located on the face of the handpiece when air passes through the generator, thus eliminating the need for a chairline fiber optic system.

- 1 Year Warranty (including generator)
- 22 Watts of Power
- Turbine With Ceramic Bearings
- Made in USA

- Surgical model provides no air into the oral cavity - all air exhausts via rear pilot holes

Built in LED Generator

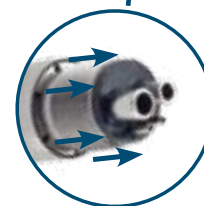

### DIRECT LED LIGHT

With bright LED optics pointed directly into your procedure site, the Foresight SGL will illuminate the procedure site better than any headlamp or loupe on the market.



### BUILT IN MICRO GENERATOR

The micro generator incorporated inside the handpiece activates a brilliant LED light when air flows through the handpiece. You no longer have to deal with inconvenient fiber optic systems.



### TITANIUM COATED STAINLESS STEEL

This titanium coated stainless steel shell provides scratch, stain, and corrosion resistance for durability and longevity.

• 22 Watts of Power - for consistent power and controlled speed for smooth rapid cutting even under heavy load conditions. • 1 Year Warranty - on the handpiece, turbine, and generator.

